




ARROW FLEXBLOCK™


Continuous peripheral nerve block catheter



TELEFLEX'S ARROW FLEXBLOCK CONTINUOUS PERIPHERAL NERVE BLOCK CATHETER

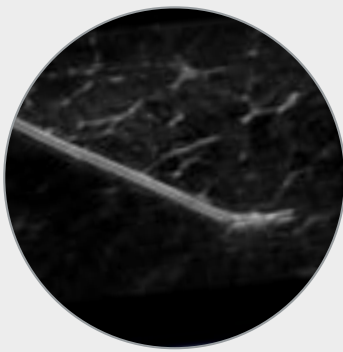

The FlexBlock catheter's coil-reinforced polyurethane body is visible under ultrasound. Image of FlexBlock catheter in porcine model using SonoSite® S-Nerve™.

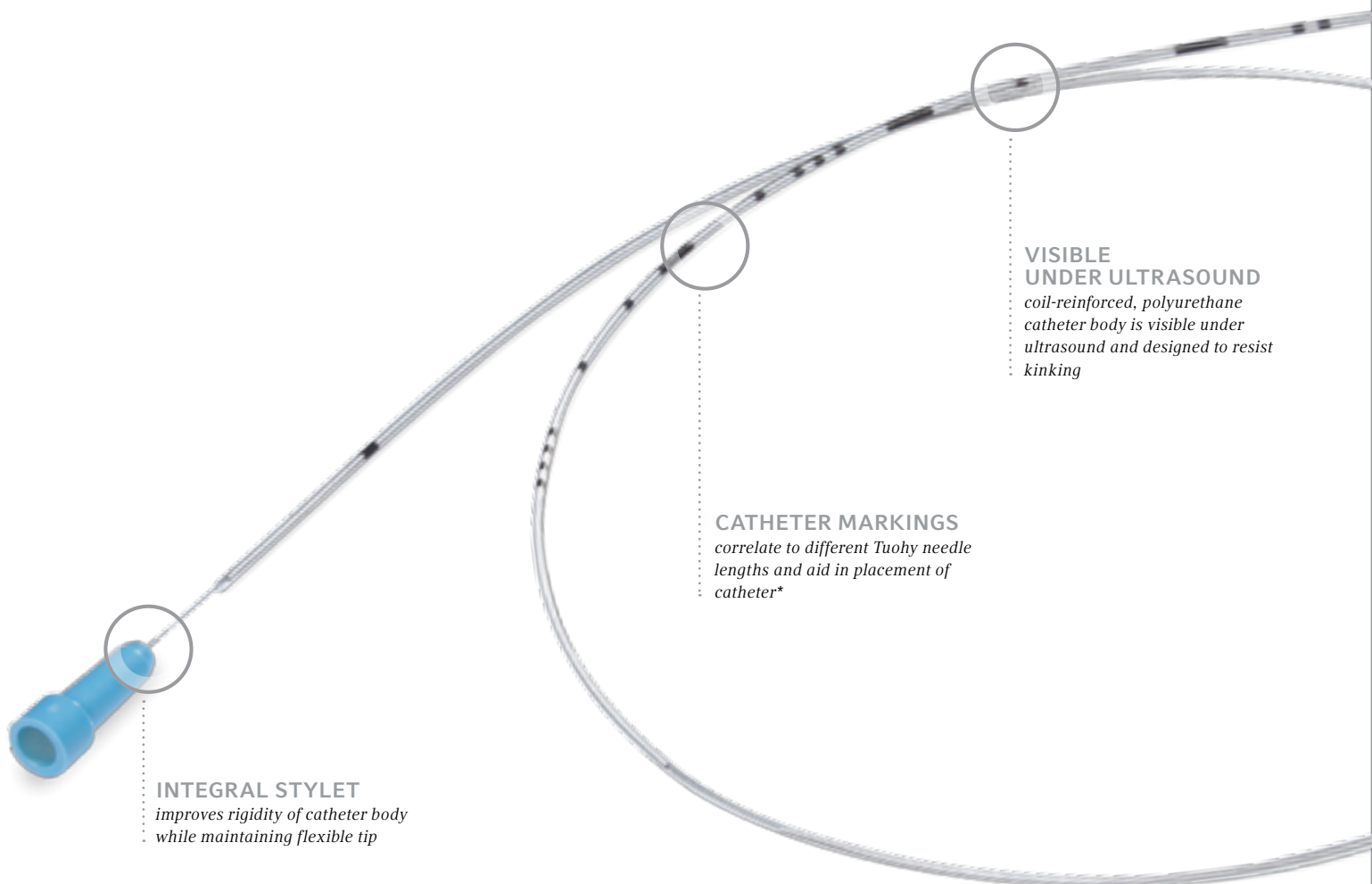


INTEGRAL STYLET

improves rigidity of catheter body while maintaining flexible tip

CATHETER MARKINGS

*correlate to different Tuohy needle lengths and aid in placement of catheter**

VISIBLE UNDER ULTRASOUND
coil-reinforced, polyurethane catheter body is visible under ultrasound and designed to resist kinking

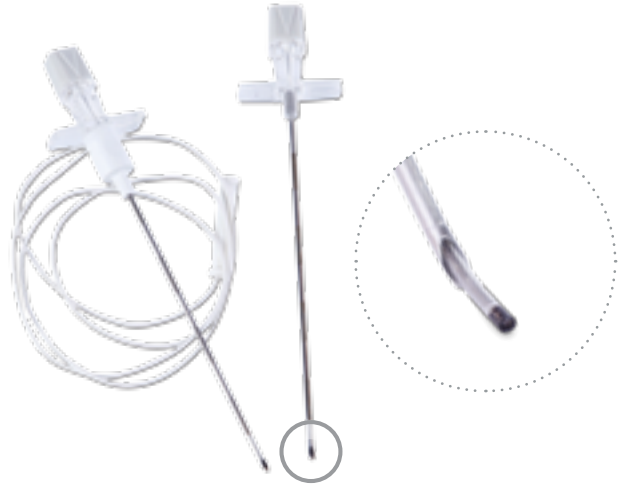


PROXIMAL CATHETER MARKING

indicates when FlexBlock catheter is properly inserted into SnapLock adapter

SPECIALLY DESIGNED TUOHY NEEDLES

- custom Tuohy needle allows for pull-back and redirection of the catheter without the need to withdraw the needle and restick the patient*
- specially-treated Tuohy needle helps increase visibility under ultrasound
- available with stimulating or non-stimulating needles

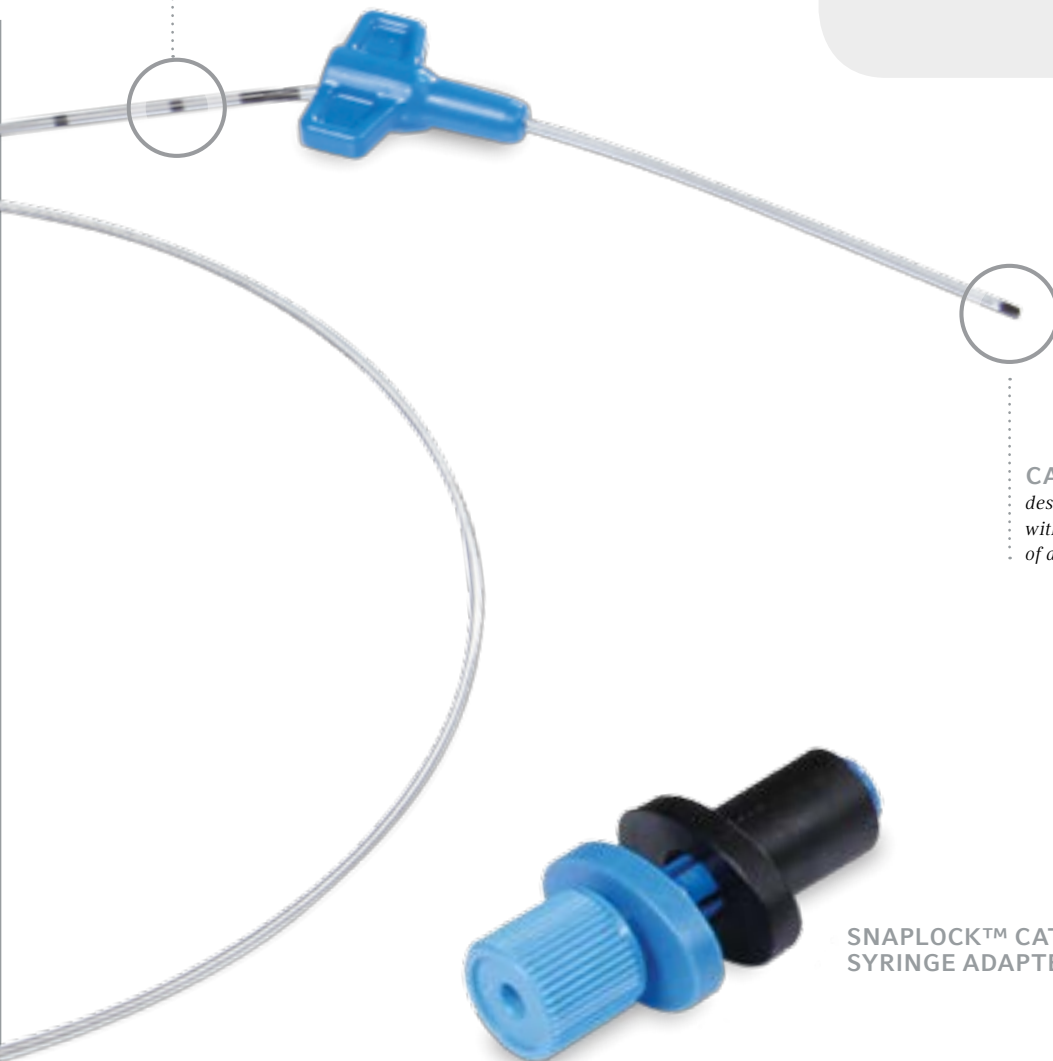

CATHETER TIP

designed to provide clinicians with a predictable spread of anaesthetic

**SNAPLOCK™ CATHETER/
SYRINGE ADAPTER**

DESIGNED TO IMPROVE OVERALL BLOCK SUCCESS

VISIBILITY

Our FlexBlock catheter's coil-reinforced, polyurethane body is designed to help improve continuous nerve block success by providing clinicians with an echogenic catheter which also resists kinking. The distal end of the catheter features a flexible catheter tip, while the integral stylet provides additional rigidity of the catheter body.

PATENCY

The coil reinforcement makes FlexBlock kink-resistant to better maintain patency, even in some situations that might cause ordinary nylon catheters to fail.

EFFICIENCY


The specially designed Tuohy introducer needles are designed to allow clinicians to retract the catheter back into the needle, if necessary, without withdrawing the needle and resticking the patient (see IFU for complete instructions).

SECURITY


All standard FlexBlock catheterisation sets include the customised StatLock catheter securement device which securely anchors the FlexBlock catheter, SnapLock adapter, and injection tubing to the patient.

THE RESULTS ARE A CONTINUOUS
PERIPHERAL NERVE BLOCK SYSTEM
DESIGNED TO HELP YOU BETTER CONTROL
YOUR PATIENTS' PAIN AND IMPROVE
OVERALL BLOCK SUCCESS.



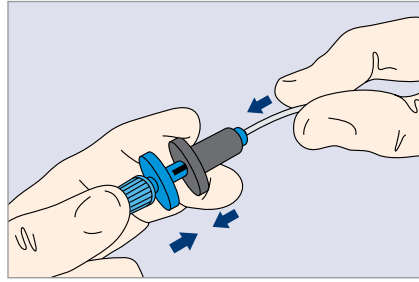

How to close the SnapLock adapter

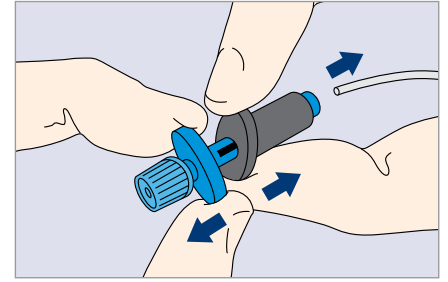


How to turn and open the SnapLock adapter

STATLOCK® CATHETER SECUREMENT DEVICE

- Custom designed StatLock anchoring device is designed to securely anchor the FlexBlock catheter, SnapLock adapter, and injection tubing to the patient

ARROW'S PROPRIETARY SNAPLOCK ADAPTER

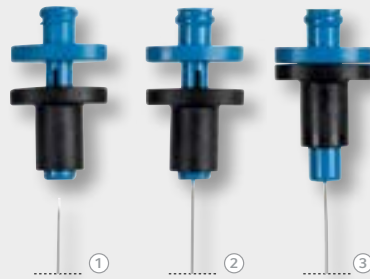

The adapter is fastened securely to the catheter, decreasing the risk of catheter disconnection.*

- ① Detached position
- ② SnapLock adapter is closed after catheter is fully inserted
- ③ Closed position

FLEXBLOCK SETS WITH STIMULATING NEEDLE

ARROW

REF.	PRODUCT	PRODUCT DESCRIPTION	QTY
FB-19604-SST	Continuous peripheral nerve block catheter with stimulating needle	<ul style="list-style-type: none"> • Tuohy needle 17 Ga. x 4 cm • catheter 19 Ga. 60 cm • Snaplock adapter, StatLock fixation tool 	10
FB-19608-SST		<ul style="list-style-type: none"> • Tuohy needle 17 Ga. x 8 cm • catheter 19 Ga. 60 cm • Snaplock adapter, StatLock fixation tool 	10
FB-19611-SST		<ul style="list-style-type: none"> • Tuohy needle 17 Ga. x 11 cm • catheter 19 Ga. 60 cm • Snaplock adapter, StatLock fixation tool 	10

FLEXBLOCK SETS WITH NON-STIMULATING NEEDLE

ARROW

REF.	PRODUCT	PRODUCT DESCRIPTION	QTY
FB-19605-SNS	Continuous peripheral nerve block catheter with non-stimulating needle	<ul style="list-style-type: none"> • Tuohy needle 17 x 5 • catheter 19 Ga. 60 cm • Snaplock adapter, StatLock fixation tool 	10
FB-19609-SNS		<ul style="list-style-type: none"> • Tuohy needle 17 x 9 • catheter 19 Ga. 60 cm • Snaplock adapter, StatLock fixation tool 	10
FB-19612-SNS		<ul style="list-style-type: none"> • Tuohy needle 17 x 12 • catheter 19 Ga. 60 cm • Snaplock adapter, StatLock fixation tool 	10

SonoSite and S-Nerve are trademarks or registered trademarks of SonoSite, Inc.
 StatLock is a registered trademark of C. R. Bard.
 Hi-Lite Orange and Chloraprep are registered trademarks of CareFusion.
 Tegaderm is a registered trademark of 3M Company.
 Teleflex, Arrow, FlexBlock, SnapLock, SnapTear, and SharpsAway are trademarks or registered trademarks of Teleflex Incorporated or its affiliates.

* See IFU for complete instructions

Teleflex is a leading global provider of specialty medical devices used for diagnostic and therapeutic procedures in critical care, urology and surgery. Our mission is to provide solutions that enable healthcare providers to improve outcomes and enhance patient and provider safety. We specialise in devices for general and regional anaesthesia, cardiac care, respiratory care, urology, vascular access and surgery and we serve healthcare providers in more than 140 countries. Teleflex also provides specialty products for medical device manufacturers.

Our well known brands include **ARROW**®, **DEKNATEL**®, **GIBECK**®, **HUDSON RCI**®, **KMEDIC**®, **LMA**™, **PILLING**®, **PLEUR-EVAC**®, **RÜSCH**®, **SHERIDAN**®, **TAUT**®, **TFX OEM**®, **VASONOVA**™ and **WECK**®, all of which are trademarks or registered trademarks of Teleflex Incorporated.

Teleflex global operations: Australia, Austria, Belgium, Canada, China, Czech Republic, France, Germany, Greece, India, Ireland, Italy, Japan, Malaysia, Mexico, Netherlands, New Zealand, Portugal, Singapore, Slovak Republic, South Africa, Spain, Switzerland, United Kingdom, Uruguay and USA.

YOUR INTERNATIONAL CONTACTS:

TELEFLEX HEADQUARTER INTERNATIONAL, IRELAND

Teleflex Medical Europe Ltd., IDA Business and Technology Park,
Dublin Road, Athlone, Co Westmeath
Phone +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73
orders.intl@teleflex.com

AUSTRALIA/NEW ZEALAND +61 (0)3 84 12 03 21

AUSTRIA +43 (0)1 402 47 72

BELGIUM +32 (0)2 333 24 60

CHINA (SHANGHAI) +86 (0)21 6381 6399

CHINA (BEIJING) +86 (0)10 6418 5699

CZECH REPUBLIC +420 (0)495 759 111

FRANCE +33 (0)5 62 18 79 40

GERMANY +49 (0)7151 406 0

GREECE +30 210 67 77 717

INDIA +91 (0)44-2836 5040

ITALY +39 0362 58 911

JAPAN +81 (0)3 3379 1511

NETHERLANDS +31 (0)88 00 215 00

PORTUGAL +351 22 541 90 85

SINGAPORE +65 6447 5022

SLOVAK REPUBLIC +421 (0)3377 254 28

SOUTH AFRICA +27 (0)11 807 4887

SPAIN +34 918 300 451

SWITZERLAND +41 (0)31 818 40 90

UNITED KINGDOM +44 (0)1494 53 27 61

For detailed information see www.teleflex.com

The products in this catalogue may not be available in all countries.
Please contact your local representative. All data current at time of printing (07/2013).
Subject to technical changes without further notice.

94 11 13 - 00 00 01 · REV A · AH / GH · 07 13 01