



## ARROW ULTRAQUIK

Peripheral nerve block needles

### AT A GLANCE

- five grooved rings at the distal end of the needle cannula to help clinicians identify the tip of the needle under ultrasound
- uninsulated needle shaft for ultrasound-only needle guidance
- atraumatic B-bevel
- easy-to-read centimetre markings
- ergonomic needle hub
- extended injection tubing length
- available in four lengths:
  - 50 mm
  - 70 mm
  - 100 mm
  - 150 mm

### CLINICIAN-INSPIRED NEEDLE DESIGN



ARROW UltraQuik peripheral nerve block needles are designed to improve the overall block success of ultrasound-guided procedures. The clinician-inspired UltraQuik needles are engineered to raise the accuracy of single-shot peripheral anaesthesia delivery by enhancing visibility under ultrasound and control and manoeuvrability during needle positioning.

To make it as straightforward as possible for clinicians to identify and guide the needle under ultrasound, the needle shaft is uninsulated and the distal end of the needle cannula has five grooved rings. There are also centimetre markings on the shaft to help confirm insertion depth visually. ARROW nerve block needles are further

distinguished by their ease of use, with an atraumatic B-bevel providing optimal balance between tactile feel and penetration, and the ergonomic needle hub being particularly easy to hold. Last but not least, the extended insertion tubing length helps reduce the risk of contamination.



## ENHANCED FEATURES OF THE ARROW ULTRAQUIK PNB NEEDLE

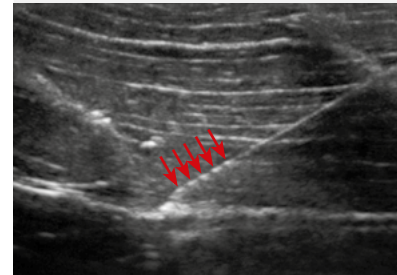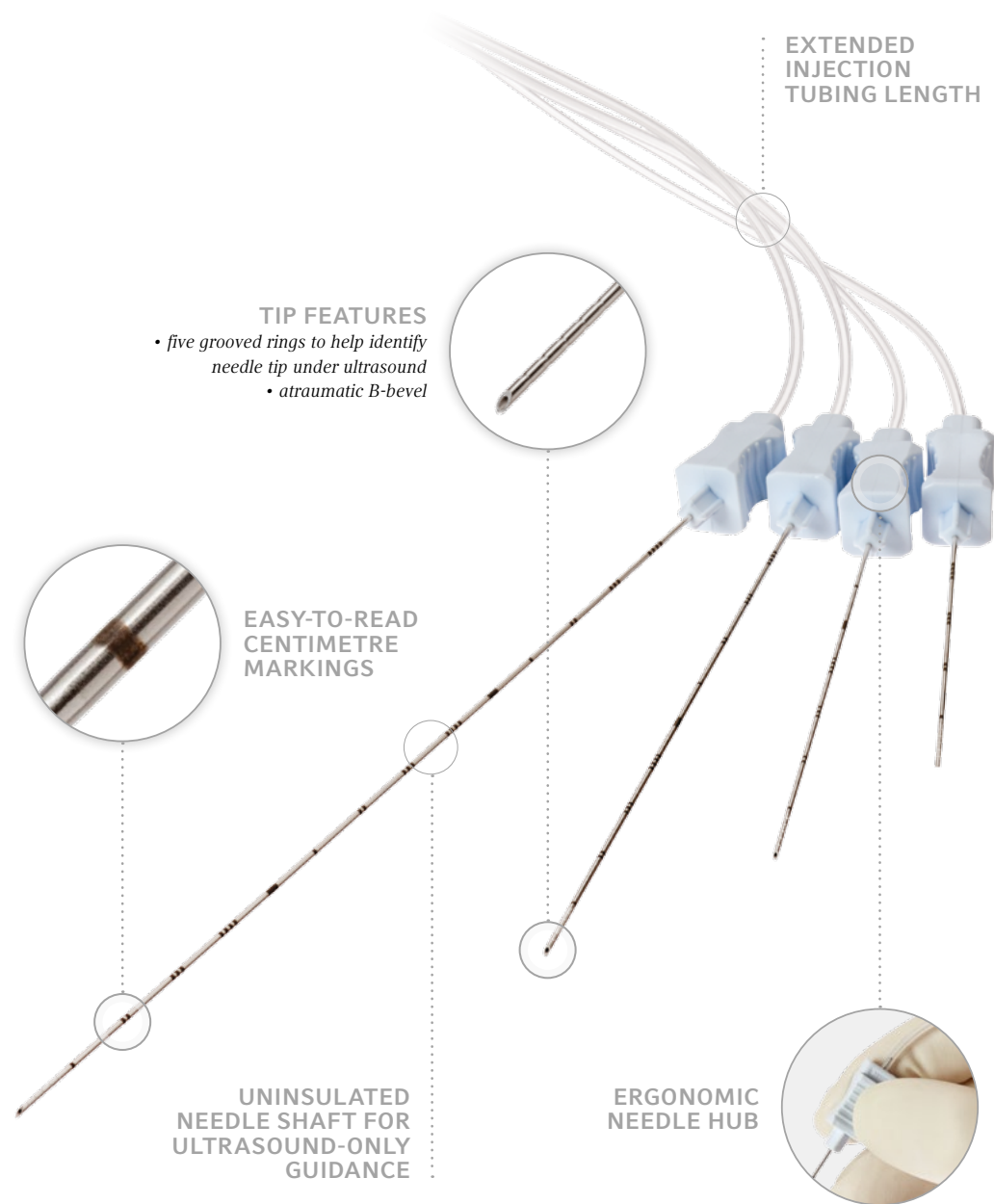

UltraQuik peripheral nerve block needle under ultrasound in porcine model using SonoSite® S-Nerve™ machine



UltraQuik features five grooved rings to help identify needle tip under ultrasound

ULTRAQUIK PNB NEEDLE			ARROW
REF.	PRODUCT DESCRIPTION		QTY
AB-22050-US	22 Ga.	50 mm	25
AB-22070-US	22 Ga.	70 mm	25
AB-21100-US	21 Ga.	100 mm	25
AB-21150-US	21 Ga.	150 mm	25

**DISTRIBUTED BY:**  
**TELEFLEX HEADQUARTER INTERNATIONAL, IRELAND** · Teleflex Medical Europe Ltd. · IDA Business & Technology Park  
 Dublin Road · Athlone · Co Westmeath · Tel. +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73 · orders.intl@teleflex.com  
**UNITED KINGDOM** Tel. +44 (0)14 94 53 27 61 · info.uk@teleflex.com  
**WWW.TELEFLEX.COM**  
 94 11 18 · 00 00 01 · REV A · 06 13 02 · All data current at time of printing. Subject to technical changes without further notice.

SonoSite and S-Nerve are trademarks or registered trademarks of Sonosite, Inc. Teleflex, Arrow and UltraQuik are trademarks or registered trademarks of Teleflex Incorporated or its affiliates.

©2013 Teleflex Incorporated. All rights reserved.